




ESSE 3 \$\$\$

Blood/Fluid warmer

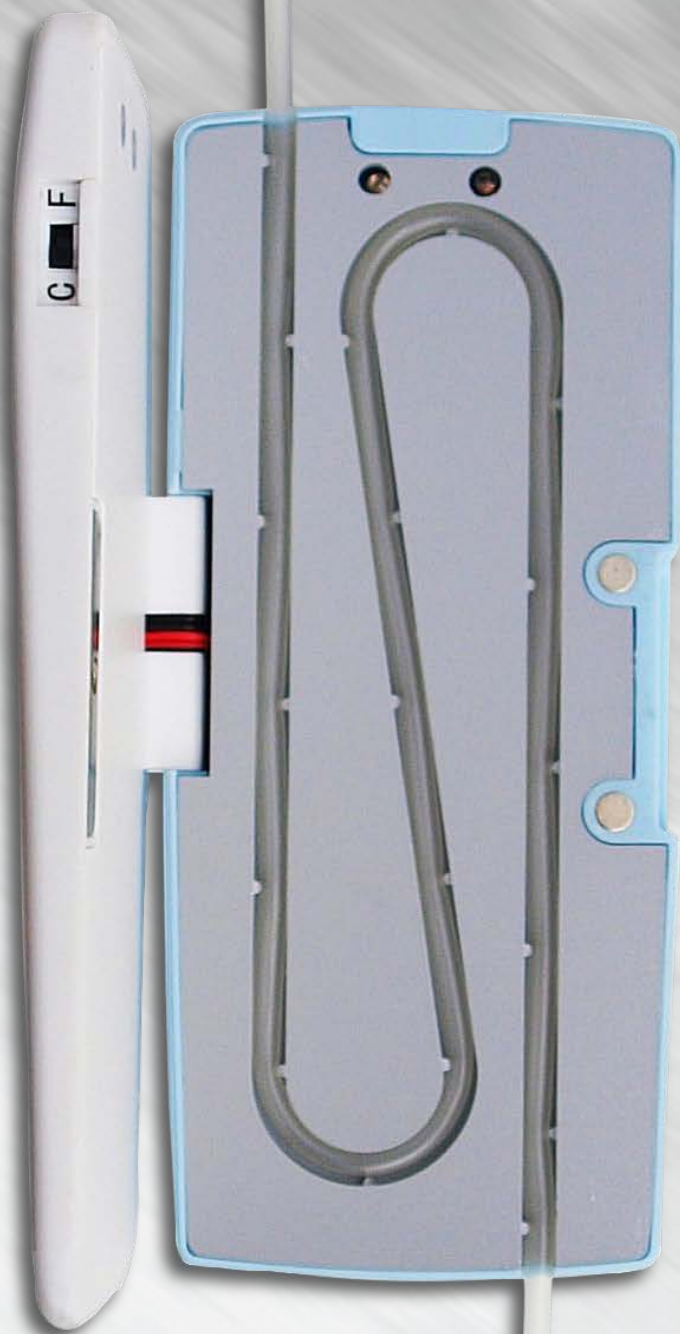

BLOOD BANK



ES<-300

...6`ccX#: `i jX`k Ufa Yf



- Completely portable & lightweight,
- Upgraded hanging clamp.
- Real-time temperature displayed in F or °C
- Adjustable temperature setting (High or Medium) for different flow rate
- Dry heat exchange warming technology, warm up any standard IV set.

Benefits for patients

- Relieves pain and disc comfort associated with the infusion of cold Fluids
- Prevents the risk of hyperthermia, peripheral vasoconstriction, cardiac arrhythmias and arrest and bleeding complications
- Prevents nosocomial infections from waterborne pathogens
Note: water is not used to warm up fluids.

Benefits for clinics/hospitals

- No need for special tubing, compatible with standard IV sets
- Convenient design and simple to use:
- Light weight and completely portable
- Powered by an A ESB-1000 infusion pump or used separately with an adaptor
- Hanging clamp included



Technical Features

Compatible IV sets	Standard IV sets
Temperature Set	Medium, High
Temperature of Heating Plate	42°C ± 5°C (Medium), 52°C ± 5°C (High) Actual fluid temperature depends on flow rate
Power Consumption	12VDC, 1.5A
Power Supply	AC 100-240V 50/60 Hz
Other Power Source	Can be powered by ESB-1000 infusion pump
Dimensions	150 x 65 x 33 mm (WxDxH)
Weight	150g / 0.33lb



ESSE3 srl, Via Garibaldi 30
14022 Castelnuovo D.B. (AT)
Tel +39 011 99 27 706
Fax +39 011 99 27 506
e-mail esse3@chierinet.it
web: www.esse3-medical.com

