


## FLOWMETERS for OXYGENTHERAPY

Mod. **PRS**

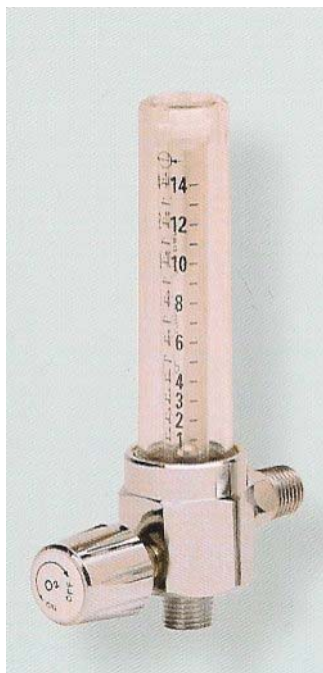


Single

The flowmeters type RS are instant flow measuring devices working at a compensated pressure particularly suitable for the use in the medical field. Complete with an adjusting needle valve, has a chrome-plated brass body, The measure group is made in polycarbonate.

### *Technical features*

- **Height:** 136 mm
- **Width:** 33 mm
- **Depth:** 84 mm
- **Weight:** Kg 0,270
- **Fluid:** O<sub>2</sub> - AIR - N<sub>2</sub> - CO<sub>2</sub>
- **End of scale value - l/min. - 3.5 bar:**  
5 - 10 - 15 - 20 - 30
- **Inlet connections:**  
ISO G 1/4"M.  
1/4" NPT M.  
3/8" ISO 3253 F.  
M12x1 F.
- **Outlet connections:**  
1/4" ISO 3253 M.  
3/8" ISO 3253 M.  
M12x1,25 M.  
1/2" GAS M.  
9/16" - 18 UNF M.  
1/2" BSF F.

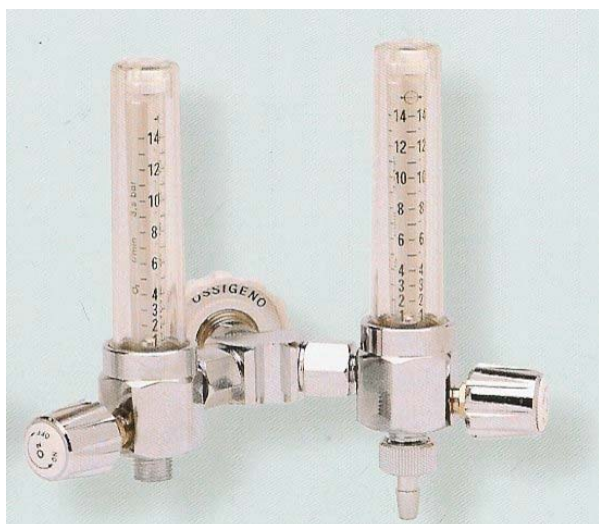

Doppio



COSTRUZIONI MEDICALI E TERAPEUTICHE


### *Laboratorio e Assistenza:*

Via Garibaldi 30 - 14022 castelnuov DB (AT)  
Tel: 039 - 011 99 27 706 Fax: 039 - 011 99 27 506  
www.4ci.com E-mail info@4ci.com



## HUMIDIFIERS FOR OXYGENTHERAPY

Mod. **TR/200** - **MAK/300** - **MAK/500**

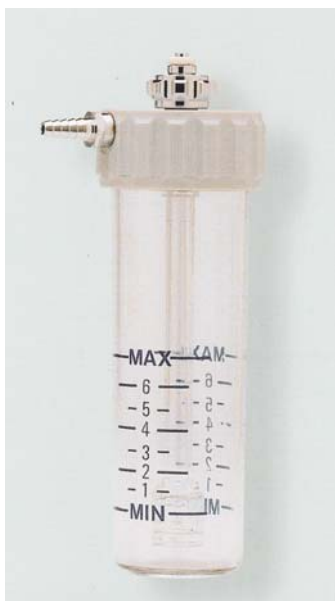

TR/200

Height 115 mm

Inner volume 200 ml



MAK/500

Height 210 mm

Inner volume 500 ml



MAK/300

Height 175 mm

Inner volume 300 ml

The bubbling humidifiers for oxygentherapy, are devices allowing to increase the relative humidity in the oxygen supplied to the patient.

Made of polycarbonate and with chrome-plated brass nipples, the oxygentherapy humidifiers are autoclavable at 134 °CC and they grant, together with the flowmeters type PRS, a simplicity in use and a versatility linked to extreme economy



**ESSE 3**

COSTRUZIONI MEDICALI TERAPEUTICHE

### Laboratorio e Assistenza:

Via Garibaldi 30 - 14022 Castelnuovo DB (AT)  
Tel: 0039 - 011 99 27 706 Fax: 0039 - 011 99 27 506  
www.4ci.com E-mail info@4ci.com

