

ES-PCA



PCA Syringe Pump


P.C.A. Function

Each PCA (Patient-Controlled Analgesia) dosage is 1/4 of preset delivery dosage.

Convenient AA Battery Power Source

Use DC 4.5V (3 AA batteries)



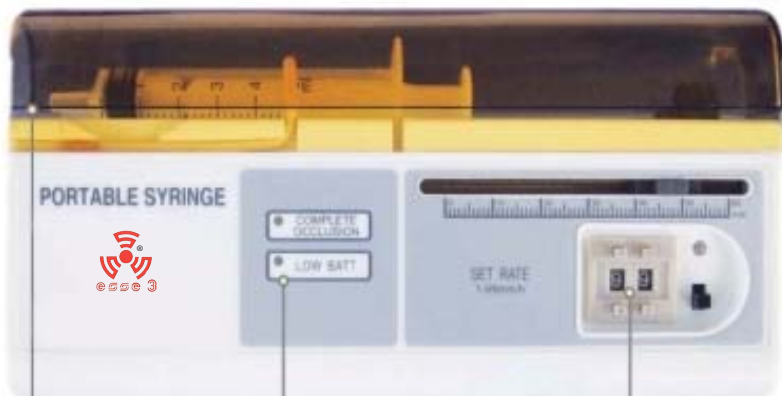

Portable/light in weight



Convenient carrying leather case & PCA bolus cable

Convenient linear-rate setting
1 to 99 mm/hr

Audible and Visual Alarms
Audible and Visual Alarms of "COMPLETE / OCCLUSION" and "LOW BATT"

Syringe protecting cover for enhanced safety

| | |
|-----------------------------|--|
| Applicable Syringes | 2ml - 10ml disposable syringes |
| Flow Rate | 1 mm/hr - 99 mm/hr (equals to 0.2 ml/hr - 20 ml/hr) |
| Syringe Accuracy | less than $\pm 2\%$ |
| Alarms | Audible and Visual Alarms of "COMPLETE / OCCLUSION" and "LOW BATT" |
| P.C.A. Function | Each PCA (Patient-Controlled Analgesia) dosage is 1/4 of preset delivery dosage. The interval between two P.C.A injections is no less than 30 minutes |
| Included Accessories | a) 4 AA batteries b) PCA bolus cable c) Syringe protecting cover d) Carrying belt |
| Dimensions | D x W x H: 30 x 166 x 60 mm D x W x H: 1.18" x 6.54" x 2.36" |
| Weight | 250 g / 0.55 lb |



ESSE 3 Via Garibaldi 30
14022Castelnuovo D.B. (AT)
tel +39 011 99 27 706
fax +39 011 99 27 506
e-mail esse3@chierinet.it
web : www.esse3-medical.com

