




ITALY & GERMANY

PRODUCT DESCRIPTION

P 750

SUTION & IRRIGATION PUMP

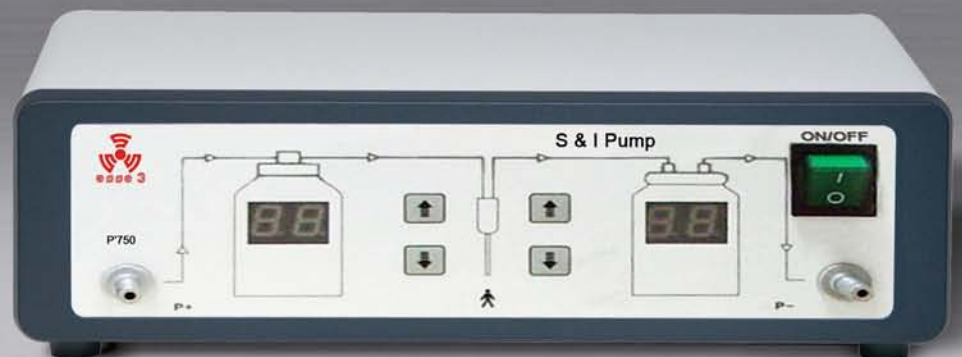

[www.4ci.com](http://www.4ci.com)

## Suction & Irrigation pump P750

The pump is a powerful suction & irrigation pump for all surgical applications and is popular using in medical operation. The system works fully automatically. During the work, no adjustments and necessary. Overpressure in the irrigation bottle and underpressure in the suction bottle are controlled and balanced by high sensitive electronical sensors and microcontroller. the automatic stand-by function keeps the pressures available and reduces the noise if the system is not used.

The use of disposable sterile air filters avoids contamination of the pump and irrigation liquid by contaminated air.

All accessories like tubes ,suction bottle and puncture tube are full autoclave.

The pump is very simple to install and easy to use.

## TECHNICAL DATA

Power supply:	AC220±10% switchable
Irrigation:	max: 3000ml/min (1)
Suction:	max: 5000ml/min (1)
Dimension :	410X415X230MM
Weight:	7kg
Working noise:	≅ 55DB
Class/type :	I/BF
Working atmosphere:	0C°-40 C°
Manufactured and tested according: ISO13485&ISO9000,CE	

## Delivery contains

P750 pump

2 liter suction bottle with holder

2 filter

1 silicone tube set for suction

2 silicone tube set for irrigation

1 puncture tube for rinsing bottle

1 Electrical source line



**ESSE 3 Export SRL**  
Via Garibaldi 30  
14022 castelnuovo DB  
Tel +39 011 99 27 706  
Fax +39 011 99 27 506  
[www.4ci.com](http://www.4ci.com)