


PRODUCT DESCRIPTION



ITALY & GERMANY

# XD 300

300 WATTS XENON COLD LIGHT SOURCE

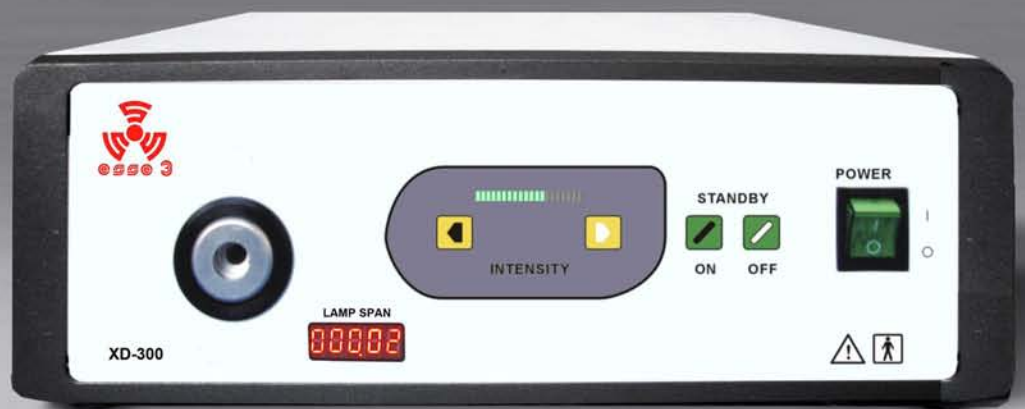

[www.4ci.com](http://www.4ci.com)

## XENON XD 300.

XD 300 offers the professional user the most advanced pure day light quality sources(6000°K color temperature). It contains an improved high intensity and low heat conductive xenon lamp.

### Lamp:

The output intensity is controlled stepless via a motor driven mechanical diaphragm. It ensures a constant color temperature. The output intensity is set by push buttons.

### Connection:

An adapter for ESSE 3 & Wolf. rigid endoscopes application.

### Quality:

The light source fulfills all essential safety demands. The rugged desk housing and the flat, waterproof touch control panel withstands all daily used requirements.

### Brightness indication:

With the LED indicator for the brightness, makes brightness more accurate.

### Lamp Span indication:

With LED lamp span indicator, the Xenon Lamp life can be easily recognized.

## TECHNICAL DATA

Model: XD-300.

Power supply: AC(110-220V)±10%,50-60HZ

Dimensions: 370x49x145

Weight: 12KG

Power consumption: 250WATT

Light guide adapter: ESSE 3 & Wolf

Life span: 500 hours

Brightness: 60000001x

Lamp: Xenon light lamp

Manufactured and tested: ISO9001/2000, CE



**ESSE 3 Export SRL**  
Via Garibaldi 30  
14022 castelnuovo DB  
Tel +39 011 99 27 706  
Fax +39 011 99 27 506  
[www.4ci.com](http://www.4ci.com)