




ITALY & GERMANY

# XD300+HL

300Watt Xenon light source with  
250Watt Halogen

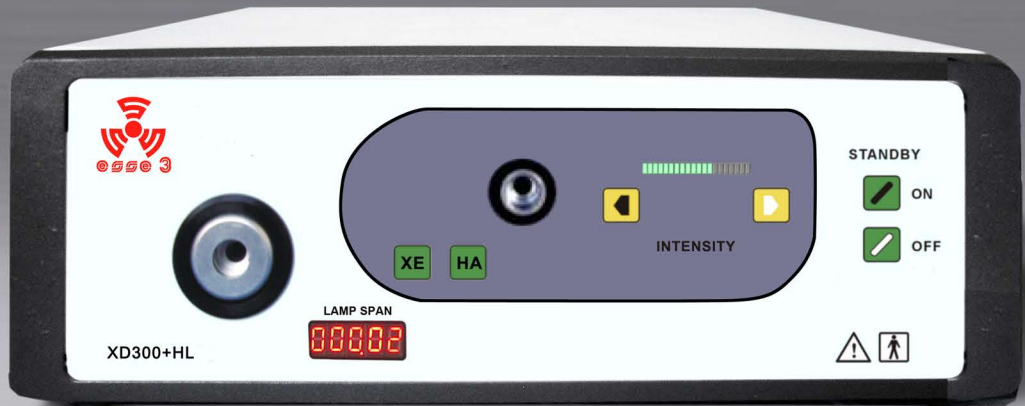

## **300Watt Xenon light source with 250Watt halogen**

XD 300 offers the professional user the most advanced pure day light quality sources(6000°K color temperature). It contains an improved high intensity and low heat conductive xenon lamp. With switching button to select 300Watt Xenon and 250Watt Halogen .

### **Lamp:**

The output intensity is controlled stepless via a motor driven mechanical diaphragm. It ensures a constant color temperature. The output intensity is set by push buttons.

### **Connection:**

An adapter for ESSE 3 & Wolf rigid endoscopes application.

### **Quality:**

The light source fulfills all essential safety demands. The rugged desk housing and the flat, waterproof touch control panel withstands all daily used requirements.

### **Brightness indication:**

With the LED indicator for the brightness, makes brightness more accurate.

### **Lamp Span indication:**

With LED lamp span indicator, the Xenon Lamp life can be easily recognized.

## **TECHNICAL DATA**

**Model:** XD300+HL.

**Power supply:** AC(110-220V)±10%,50-60HZ

**Dimensions:** 370x49x145

**Weight:** 12KG

**Power consumption:** 250WATT

**Light guide adapter:** ESSE 3 & Wolf

**Life span:** 500 hours

**Brightness:** 60000001x

**Lamp:** 300Watt Xenon & 250Watt Halogen

**Manufactured and tested:** ISO9001/2000, CE

**Light guide selected :**

**900209 : 2.5meters fluid light guide**

**900210 : 3 meters fluid light guide**

**ESSE 3 Medical**  
Via Garibaldi 30  
14022 castelnuovo DB  
Tel +39 011 99 27 706  
Fax +39 011 99 27 506  
[www.esse3-medical.com](http://www.esse3-medical.com)